


# Bed Alarm – Integrated, Stryker MED-Surg Gobed

	<p><b>Indications:</b></p> <p>The Bed Exit system is a falls prevention strategy intended only to assist in the detection of a patient exiting the bed.</p>
	<p><b>Considerations:</b></p> <p>The bed exit system is not designed to be used with patients weighing less than 50 lb (23 kg)</p> <p>The alarm signal will be heard only in the patient's room</p> <p>The Bed Exit system <b>MUST</b> be zeroed before use.</p> <p>The addition or removal of equipment with the bed exit system armed must be done using the "Adding or Removing Equipment with the System Armed" procedure, otherwise the sensitivity of the system may be affected and the readings of the patient's movement in the bed be erroneous.</p>
 <p><b>Figure 2.14C</b></p> <p>Press Arm/disarm (K) to activate or deactivate the Bed Exit</p> <p>LED status light will display (M) when armed</p>	<p><b>Arming the System</b></p> <p>Before putting a new patient on the bed, prepare the bed for a patient stay by adding linens and equipment to the bed. The Bed Exit system absolutely must be zeroed before use. Following the zeroing of the system, the Bed Exit system may be activated.</p> <ul style="list-style-type: none"> <li>Position the patient in the bed. Ensure the patient lies in the bed in a usual position before activating the system.</li> <li>Activate the Bed Exit by pressing Arm/Disarm (K). The status LED (M) will light up.</li> <li>To deactivate a Bed Exit alarm or to turn the system off, press Arm/Disarm (K). The status LED (M) will turn off. The Bed Exit system may be activated once the patient returns to the bed.</li> </ul>

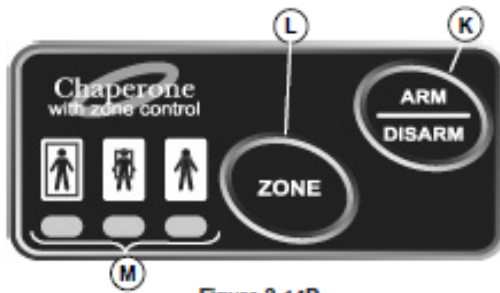


Figure 2.14D

Press Arm/Disarm (K) to activate/deactivate the Bed Exit

Press Zone to select the required monitoring zone

The LED status light will indicate the zone selected (M)

MODEL: FL 28C

#### BED EXIT SYSTEM WITH MULTI-ZONE CONTROL

Failure to set the zone properly could result in improper monitoring and patient injury. Verify that the proper zone is selected before leaving the patient.

##### Zone Features

- The first zone (left LED) is the traditional Bed Exit zone. The patient can move in the bed freely but cannot fully exit the bed or the alarm will sound.
- The second zone (middle LED) is more restrictive. It allows the patient to sit up and roll over but any attempt to exit the bed will cause the alarm to sound.
- The third zone (right LED) is the most restrictive. Small movements like raising an arm or lifting the shoulder off the bed will cause the alarm to sound. The third zone is used to alert staff to a change in the condition of an unconscious or paralyzed patient.
- Press and hold the ARM/DISARM key (K) until the larger zone LED (left LED) begins to flash, then release the control (do not touch the bed while the LED flashes). When the zone LED turns off, the system is zeroed.
- By default, the least restrictive zone will be selected and its LED (left LED) will come on. To change the zone, press ZONE (L) repeatedly until the desired zone LED comes on.
- To deactivate the Bed Exit after an alarm or to simply turn it off, press ARM/DISARM (K). The zone LED will turn off.